



## The Monaghan Advantage

- ❖ Breath actuated technology nebulizer
- ❖ Durable, 6-month reusable design
- ❖ Allows for breaks in treatment without drug loss
- ❖ Minimal exposure of caregivers and family members to secondhand aerosols
- ❖ Optimal performance with the OMBRA® Compressor System
- ❖ Nebulizer designed for high level disinfection and dishwasher safe
- ❖ Not made or manufactured with BPA, natural rubber latex, DEHP or lead

## AEROECLIPSE® XL BAN AND OMBRA® COMPRESSOR

The AEROECLIPSE® XL BAN provides a hospital-quality nebulizer treatment for the home care environment. Recommended by clinicians, the AEROECLIPSE® XL BAN is designed to improve patient compliance for nebulizer treatments outside of the hospital.

Activating only in response to the patient's inspiratory flow, AEROECLIPSE® XL BAN maintains available drug during exhalation phase, providing patients the ability to take a break or have a conversation without drug loss. Paired with the OMBRA® Compressor, the AEROECLIPSE® XL BAN puts patients in control of their therapy. Unique and patient friendly features make the system an easy and safe to use, cost effective drug delivery system.

## THE FACTS

**Recommended by clinicians, the AeroEclipse® XL BAN is designed for nebulizer treatments to improve patient compliance outside of the hospital.**

**Suitable for use for by all patient populations, the AeroEclipse® XL BAN activates at 15 lpm and can be used with a mouthpiece or our ComfortSeal® Masks.**

**Minimizes exposure of caregivers and family members to dangerous secondhand aerosols, which may carry harmful viruses.**


**Designed to put patients in control of their therapy at home, AeroEclipse® XL BAN maintains available drug during exhalation phase to give patients time to take a break or have a conversation without drug loss.**

MAINTAINS  
AVAILABLE  
DRUG DURING  
EXHALATION  
PHASE



518.561.7330 phone  
800.833.9653 toll free  
518.561.5088 fax

5 Latour Avenue, Suite 1600  
Plattsburgh, New York 12901  
[monaghanmed.com](http://monaghanmed.com)

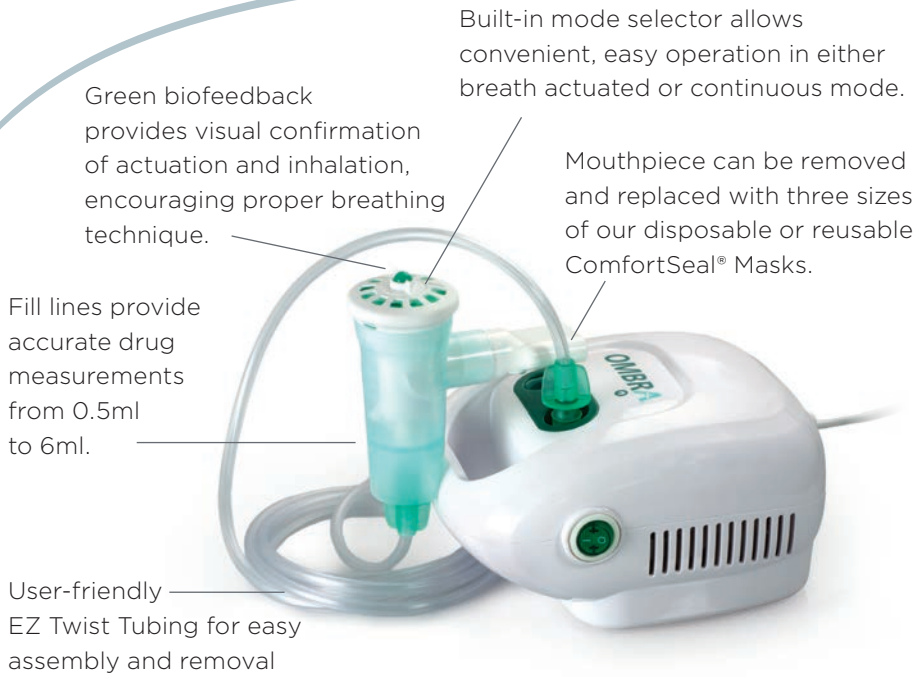

## THE FACTS (continued)

**Pair the AeroEclipse® XL BAN with the Ombra® Compressor for optimal nebulizer performance in a durable home compressor system.\* Average treatment time is 6–8 minutes.**

**Reusable for 6-months of single-patient use, the AeroEclipse® XL BAN is top shelf dishwasher safe and approved for high-level disinfection.<sup>1</sup>**

**Kit includes 2 AeroEclipse® XL BAN with integrated handle and nebulizer cup holder.**

### REFERENCE

1. TMI Aerosol Lab, Data on File. 2.5ml, 2.5 mg salbutamol sulphate (Vt = 600 mL, i:E 1:2, 10 bpm) @ 8 minutes ‡ Mass Median Aerodynamic Diameter 3.4.

\*Warranty information about our OMBRA® compressors is available upon request by contacting Monaghan Customer Service at 1-800-833-9653.

## CLEANING IS EASY

- AeroEclipse® XL BAN, single patient use
  - Designed for 6 months of use at home, clinic or hospital environment
  - Approved for top shelf dishwasher cycle and high level disinfection with bleach, boiling, microwave steam bag or electronic steam sterilizer

## ORDERING INFORMATION



PART #	NAME/DESCRIPTION	UNITS PER CASE
69010	Reusable AEROECLIPSE® XL Breath Actuated Nebulizer (BAN)	cs/10
63506	OMBRA® 120V Table Top Compressor (includes 2 AEROECLIPSE® XL BAN, 5 Replacement Air Filters, Owner's Manual)	cs/1
66501	OMBRA® Compressor Air Filters, Replacement (5/pack)	cs/1

Contact your Monaghan Sales Territory Manager for more information or a demonstration.

